



# Simplify with SimplyGo

Now one POC meets the portable oxygen needs of nearly all of your patients

**Learn More at**  
**[OxygenConcentratorStore.com](http://OxygenConcentratorStore.com)**



**OXYGENConcentratorStore**

FROM AMERICAN MEDICAL SALES & RENTALS

**1-877-774-9271 | [www.oxygenconcentratorstore.com](http://www.oxygenconcentratorstore.com)**



# Treating your oxygen patients has never been easier

SimplyGo is the only portable oxygen concentrator (POC) to offer continuous flow and pulse-dose delivery in a single device weighing only 10 pounds.

That means you can now manage the therapy and lifestyle needs of nearly all of your patients with just one POC. Designed to be rugged and easy to use, SimplyGo is truly a new breed of lightweight and attractive portable concentrator. Simplify patient management – and your business – with the SimplyGo portable oxygen concentrator from Philips Respironics.

A new breed of portable concentrator that simplifies the care and management of nearly all oxygen patients.



1-877-774-9271 | [www.oxygenconcentratorstore.com](http://www.oxygenconcentratorstore.com)



# SimplyGo simply...



## Meets most patients' needs

SimplyGo packs more of what you need in one POC, without packing on unwanted pounds. In fact, SimplyGo is the only continuous flow POC to weigh in at ten pounds. And for patients who need higher pulse settings at night, the Sleep mode provides increased trigger sensitivity and a softer pulse for exceptional comfort. There's simply no POC more flexible than SimplyGo.

## Redefines ruggedness

When it comes to toughness, SimplyGo is no pushover. With a long-life compressor, high-quality parts, a high-impact resistant design, and oversized cart wheels, this POC was built to set new standards of ruggedness. During development, SimplyGo was rigorously tested and subjected to extreme conditions including impacts, vibrations, and temperatures to deliver reliable performance day in and day out.


## Fits your business

SimplyGo's warranty and service programs are all about flexibility. With SimplyGo, you have the choice of a two-, three-, or five-year warranty. Pick the plan that's best for you. While some POCs must be returned to the manufacturer for service, SimplyGo's easily accessible parts and comprehensive service tools allow your technicians to perform the work right in your shop if you prefer.

# SimplyGo for today and tomorrow

SimplyGo has more than twice the oxygen output of any POC weighing ten pounds or less. That means that SimplyGo is able to keep up with your patients' portable oxygen needs today and tomorrow.

## Oxygen output of POCs weighing 10 pounds or less

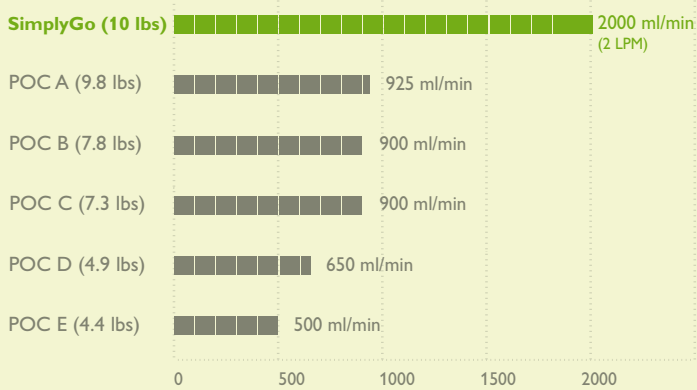

## Powered for success

One of the keys to SimplyGo's portability is its lightweight and compact battery. SimplyGo's battery is immediately accessible and can be easily changed by your patients. There are no awkward external battery belts or connectors that get in the way, or hard-to-maintain internal batteries.

### SimplyGo's battery duration

Pulse dose setting of 2*	Continuous flow 2 LPM
3.0 hours	.9 hours

\*Based on 20 breaths per minute.

Actual duration results may vary.



The basic SimplyGo system comes with the SimplyGo unit, AC/DC adapter, carrying case, shoulder strap, mobile cart, bag, cannula and manual. Battery not included. But if you want more power, extra batteries fit easily into the zippered pouch or handy accessory bag. Patients want more time? Simply take along more batteries.

## Intended to streamline your business...

Since SimplyGo meets the needs of nearly all of your patients, you can greatly simplify your inventory management and use less shelf space, helping you to lower costs and streamline patient services. And, SimplyGo's ability to provide continuous flow in a small, lightweight package is appealing to referral sources and patients.



## ...and be a team player

When SimplyGo is teamed up with EverFlo – our low maintenance, lightweight, stationary oxygen concentrator for use in the home – you can create an oxygen-generating system that will meet the needs of nearly all of your patients while helping to reduce your cost of providing oxygen therapy.

Our SimplyGo/EverFlo duo helps:

- Eliminate recurring oxygen delivery costs and cylinder hassles
- Reduce maintenance costs
- Ease inventory management
- Minimize patient calls



## One attractive package

SimplyGo's non-clinical style and ease-of-use make it appealing to patients. An attractive and highly functional carrying case allows users to either carry SimplyGo over the shoulder or place it on the easy-to-pull, easy-to-store cart. Power from a Lithium Ion battery or standard AC/DC power cords add convenience and flexibility. So, whether it's a walk across the street, a ride across town, or a trip across the country, we've made it easier than ever for you to help patients – and yourself – simplify with SimplyGo.



## The SimplyGo portable oxygen concentrator

Comes with everything that you need including the SimplyGo unit, AC/DC adapter, carrying case, shoulder strap, mobile cart, bag, cannula and manual. Battery sold separately.

- Attractive, sleek design
- Intuitive, easy-to-read screen
- Convenient carrying case
- Handy fold-up cart
- Detachable accessory bag



Airline policies for traveling with a portable oxygen concentrator vary. Patients should check with their airlines before flying with SimplyGo.

1-877-774-9271 | [www.oxygenconcentratorstore.com](http://www.oxygenconcentratorstore.com)

Specifications	
Oxygen concentrations	87 to 96% at all settings
Pulse mode bolus size	1=12 ml, 1.5=18 ml, 2=24 ml, 2.5=30 ml, 3=36 ml, 3.5=42 ml, 4=48 ml, 4.5=54 ml, 5=60 ml, 5.5=66 ml, 6=72 ml
Continuous flow settings	0.5 LPM, 1 LPM, 1.5 LPM, 2 LPM
Sound level	≤ 43 dB at pulse setting 2
Size	11.5 in x 10 in x 6 in (29.2 cm x 25.4 cm x 15.2 cm)
Weight	10 lbs (4.5 kg) with battery
Operating temperature	41° to 104°F (5° to 40°C)
Storage temperature	-4° to 140°F (-20° to 60°C) unit only
Altitude	Up to 10,000 ft. (3048 m)
Filtration	Internal bacteria filter
DC power	Input 11.5 to 16 VDC, 18 A max Output 19 VDC, 7.9 A max
AC power	Input 100 to 240 VAC, 50/60 Hz, 1.7 A max Output 19 VDC, 7.9 A max
Power consumption	150W while charging, 120W while not charging
Alarm types	Low oxygen purity, No breath, High breath rate, No flow, Low battery, Depleted battery, Wrong battery, System technical fault
Lithium Ion battery	14.4 VDC, 6600 mAh (Lilon), 1.5 lbs (0.7 kg), charge time is approximately 2 to 3 hours per battery from fully discharged

Ordering information	
<b>Basic system</b>	
Includes SimplyGo portable oxygen concentrator, carrying case, mobile cart, DC and AC power cords, accessory bag, cannula and user manual	1068987
<b>Accessories</b>	
SimplyGo carrying case	1082663
Lithium Ion battery	1082662
Battery charger	900-103
DC power supply	1083692
AC power supply and cord	1082661
Airline power cord	1083693
Accessory bag	1083696
Mobile cart	1074885
User manual	1083700
Humidifier pouch	1083699
SimplyGo humidifier connecting tube kit	1101172
SimplyGo humidifier connecting tube kit and pouch	1101602
<b>Warranty options</b>	
SimplyGo two-year total warranty	1097623
SimplyGo three-year total warranty	1097577
SimplyGo five-year total warranty	1097778



Please visit [www.oxygenconcentratorstore.com](http://www.oxygenconcentratorstore.com)



**OXYGENConcentratorStore**

FROM AMERICAN MEDICAL SALES & RENTALS

**1-877-774-9271 | [www.oxygenconcentratorstore.com](http://www.oxygenconcentratorstore.com)**