



A new choice for nocturnal patients...and your business

Introducing **SimplyFlo**, the world's smallest and lightest stationary concentrator



1-877-774-9271 | www.oxygenconcentratorstore.com

A new solution for nocturnal patients

Once again, Philips Respironics is bringing innovative solutions to your business and your patients. With SimplyFlo, we've created the world's smallest and lightest stationary oxygen concentrator.

SimplyFlo represents a completely new option for nocturnal patients who may be refusing or resisting therapy because of the size, weight and lifestyle disruptions that can accompany the use of traditional concentrators.

Designed to enhance compliance

For many patients new to therapy who only require nocturnal treatment, the intimidating size, weight, and sound of a traditional oxygen concentrator can be a large factor in their rejection of therapy.

Until SimplyFlo, there has not been a significantly lighter, more transportable stationary concentrator designed to help nocturnal oxygen patients improve the adoption of therapy and enhance compliance.



The optional travel accessory case allows patients to easily transport SimplyFlo between locations.





PHILIPS
RESIRONICS

A new opportunity for your business




In our recent survey of homecare providers, 75% reported having to pick up or accept the return of a traditional concentrator due to patient refusal of therapy.*

*Data on file.

Your nocturnal-use oxygen patients – and your business – require an option to help patients who may be refusing therapy. SimplyFlo is designed to be that option.

Compliance among nocturnal patients

Homecare providers report a range of reasons for rejection of therapy among new nocturnal patients:*

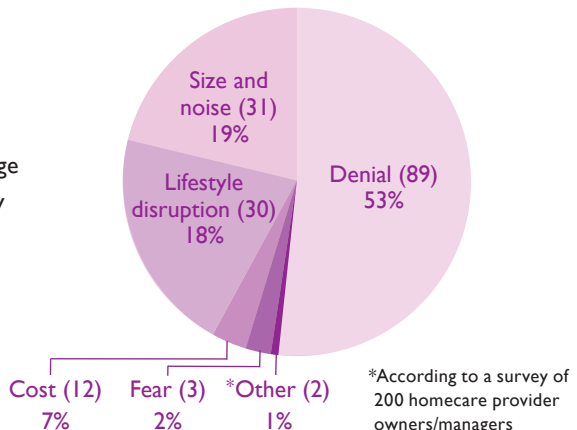

Benefits ideal for nocturnal patients



- Reduces many of the lifestyle disruptions related to larger stationary concentrators
- Simplifies transport for travel or within home
- Helps nocturnal oxygen users maintain an active lifestyle during the day

SimplyFlo features

Weights 8.5 pounds

Provides .5 to 2 LPM continuous flow output

Operates on AC power (no internal battery)

SimplyFlo is a great opportunity for new retail sales to your existing nocturnal patients as a convenient, transportable secondary device.



SimplyFlo uses less than 1/3 the electricity of traditional stationary concentrators

For more information about SimplyFlo, please contact your American Medical Sales and Repair

1-877-774-9271 | www.oxygenconcentratorstore.com

SimplyFlo stationary oxygen concentrator

Specifications

Product specifications	
Oxygen concentrations*	87 to 96% at all settings
Outlet pressure	6.4 PSIG max
Continuous flow settings	0.5 LPM, 1 LPM, 1.5 LPM, 2 LPM +/- 15% or 150 ml/min whichever is greater (3 minute running average)
Sound level	≤50 dBA (typical) at 2 LPM
Physical characteristics	
Size	11.5 in x 10 in x 6 in (29.2 cm x 25.4 cm 15.2 cm)
Weight	<9 lbs (4.1 kg)
User Interface	Push buttons, back-lit liquid crystal display (LCD)
Filtration	Internal: air inlet filter (gross particle), main-flow bacteria filter
Operating conditions	
Operational temperature	41° to 104° F (5° to 45°C)
Relative humidity	15 to 95%
Altitude	Up to 10,000 ft. (3048 m)
Storage conditions	
Temperature	-4° to 140°F (-20°F to 60°C) unit only
Relative humidity	Up to 95%, non condensing
AC power	
Input	100 to 240 VAC, 50/60 Hz, 1.7 A max
Output	19 VDC, 7.9 A max
Power consumption	120W max
Audible alarm	
Sound level	50 dB (nominal) at 39 in (1 m)
Alarm type	Low oxygen purity. No flow. System technical fault. Loss of power

*Based on atmospheric pressure of 14.7 psi (101kPa) at 70°F (21°C)

Ordering information

SimplyFlo system	
Basic system	1076773
Includes SimplyFlo oxygen concentrator, AC power cord, and user manual	
Accessories	
SimplyFlo travel case	1076830
AC power cord	1082664
AC power supply	1082661
User manual	1076809
Humidifier bottle holder	1076829
Humidifier connecting tube kit	1101172
Warranties	
SimplyFlo two-year total warranty	1097623
SimplyFlo three-year total warranty	1097577

Reimbursement information

HCPCS code	Equipment	Description
E1390	Oxygen concentrator—single delivery port	Single delivery port, capable of delivering 85% or greater oxygen concentration at the prescribed flow rate, each



Allies
in Better Sleep
and Breathing

Philips Respironics is working in harmony with the care team to produce solutions that can lead to healthier patients, healthier practices and healthier businesses.

1-877-774-9271 | www.oxygenconcentratorstore.com

Please visit www.oxygenconcentratorstore.com



American Medical Sales and Repair
7042 S Revere Parkway
Suite 450
Centennial, CO 80112

1-877-774-9271 | www.oxygenconcentratorstore.com